

Ramvac

DENTALEZ

UTILITY ROOM ESSENTIALS

ADJUSTABLE UTILITY STAND

The stand that saves utility room space and supports your Ramvac® needs

- Easy assembly with minimal tools required
- Height adjustable in 2" increments from 40" to 62"
- Load rating is a hefty 600 pounds
- Rubber mat to reduce noise
- Wide open front and sides for easy accessibility

Product Description

Part Number

Adjustable utility stand 40-62"
(36"W x 24"D x 40-62"H)

4060



Adjustable Utility Stand, shown at 62, 50 and 40 inches.

NOTES:

The Ramvac C18-2 Compressor with sound cover and RV12 Vacuum with 30 gallon tank do not fit underneath the stand.

NXT HG5 AMALGAM SEPARATORS

Market leader in the industry

- ISO 11143 certified by NSF International

Flexible design

- Adjustable for hard-to-fit spaces with right or left orientation

Compatible with wet and dry vacuum systems

- Install before pump on wet vacuum systems
- Install before tank on dry vacuum systems

NXT Hg5 is clear by design

- Allowing visual inspection of system
- No tools necessary, no daily maintenance, no decanting

Easy to change out of the collection container

- Simple to change either by dental personnel or technician



Product Description	Part Number
NXT Hg5 Amalgam Separator (1-10 users)	8707-002
NXT Hg5 MINI Amalgam Separator (1-4 users)	8707-003
NXT Hg5 High Volume Amalgam Separator (11-20 users)	8707-004

VACUUM LINE CLEANER

SlugBuster™

- Non-Foaming
- Chelators:
 - Holds Mercury in suspension for safer discharge to waste.
 - Softer water allows for better cleaning.
- 8 ph
- Wetting Agents:
 - Specialty formulated to allow the cleaner to penetrate faster.
 - Compatible with all amalgam separators.



Product Description	Part Number
Powder Starter Pack includes a 2 gallon mixing jug, 2 40oz. containers of SlugBuster powder with scoop	900143
Powder 6 Pack includes six 40oz. containers of SlugBuster powder with scoop	900146
Liquid Single Gallon includes a single gallon of SlugBuster concentrate and hand pump dispenser	900241
Liquid Starter Pack includes 2 gallon mixing jug, 1 gallon of SlugBuster concentrate and hand pump dispenser	900243
Liquid 4 Pack Case includes 4 single gallons of SlugBuster concentrate and 4 hand pump dispensers	900246

BUCK/BOOST TRANSFORMERS


- Regulates the incoming voltage to your equipment
- Increases or lowers voltage by 24 volts

Product Description	Part Number
Facility voltage range w/o transformer 115VAC 230VAC up to 40 AMP Circuit 94-98VAC 206-212VAC	64637013
Facility voltage range w/o transformer 115VAC 230VAC up to 20 AMP Circuit 94-98VAC 206-212VAC	64637006
Auxillary Vacuum Adapter	590504



VACHECK®/FLOWCHECK®

- Check the vacuum performance for your treatment room

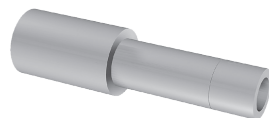

VACHECK for Simple Pass-Fail indication. (Patent 5,779,472)



FLOWCHECK for measurement of flow in SCFM. (Patent 6,925,886 B2)

Product Description	Part Number
VACHECK	800500
FLOWCHECK	800501

MASTER WATER CONTROLS

Product Description	Part Number
Master Water Control Valve, adjustable between 1/2" and 3/4"	MWC075
Master Water Control Valve 1"	MWC001
Master Water Control Valve with Regulator, adjustable between 1/2" and 3/4"	MWC075R

Replacement Cartridges	Part Number
AP110 Water Filter for 1/2" and 3/4" Master water controls (2 pack, for older WC series)	64520006
AP801 Water Filter for 1" Master water controls (2 pack, for older WC series)	64520012
Filter Cartridge Kit for new MWC series (2 pack)	8801-040

Master Water Control Valve Replacement Kits	Part Number
Solenoid Valve Kit, 3/4"	8801-041
Solenoid Valve Kit, 1"	8801-042
Transformer Kit	8801-043
O-Ring Replacement Kit	8801-044



NOTES:
Master Water Control Valves now come standard with 24V Solenoid valve and transformer for easy installation to Ramvac master control panels. 115V versions are available; please call 1-866-DTE-INFO for ordering information and availability.

FRESH AIR INTAKE



- Delivers fresh air to compressor motor heads
- Keeps compressor motor heads cool

Product Description	Part Number
Fresh Air Intake Kit, Ramvac Dual Voltage C3-12	8658-009
Fresh Air Intake Kit, Ramvac Dual Voltage C4-12	8658-010
Fresh Air Intake Kit for C6, C8, C12	8658-012
Fresh Air Intake Kit for C18	8658-013
Fresh Air Intake Kit for C3-12SS	8658-014
Fresh Air Intake Kit for C4-12SS & C5-12SS	8658-015
Fresh Air Intake Kit for C8-2SS & C12-2SS	8658-019
Fresh Air Intake Kit for C18-2SS	8658-020
Fresh Air Intake Kit for 3-7 user	005002
Fresh Air Intake Kit for 11 user	005014

REMOTE PANELS

- Control utility room equipment from front office
- Available in 1, 2, 3, 4 and 5 button configurations

Product Description	Part Number
Remote Panel #1	8-600-131
Remote Panel #2	8-600-132
Remote Panel #3	8-600-133
Remote Panel #4	8-600-134
Remote Panel #5	8-600-135

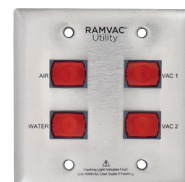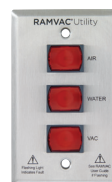

RAMVAC TANKS

- Store liquids and debris removed from the treatment room
- Tanks must be sized to correctly match the volume of liquids and debris removed during a vacuum operating cycle.

Product Description	Part Number
15 Gallon Tank**	A575015
30 Gallon Tank**	A575030

Tank Options

Ramvac Vacuum Sizes 5-12 Platform	575056
Suspension Set	575215
12" Height Extender	575212
Auxiliary Vacuum Adapter	590504

Specifications	15 Gallon	30 Gallon
Treatment Rooms*	1-6	7-12
W	29"	29"
D	19"	19"
H	24"	36"

RAMVAC INFINI TANK

- Simple maintenance free design
- Drains without vacuum interruption
- Significantly reduces the size and number of separating tanks
- Completely automated to provide worry free operation

Product Description	Part Number
Ramvac Infini Tank	A575080

Specifications	Infini Tank	Infini Tank High Flow
Treatment Rooms*	1-20	1-20
W	29"	55"
D	21"	27"
H	38"	62"

* This assumes 1 chair per room. For bays, let "rooms"= chairs.

** Ramvac tanks include level sensor and cable interface to Ramvac Vacuum Control with Aeras Intel.

NOTES:

Please refer to your User Guide or contact Ramvac Technical Support to determine your Ramvac's oil capacity.

©2023 DENTALEZ, Inc. DENTALEZ, FLOWCHECK, and VACHECK are registered trademarks of DENTALEZ, Inc. Ramvac is a registered trademark and SlugBuster is a trademark of Ramvac Dental Products Inc.

